

Size and colour chart

Arm Sleeve – Child





By Second Skin

READY TO WEAR GARMENTS

SIZES

Size	Upper arm circumference (a)	Elbow circumference (b)	Wrist circumference (c)	Length axilla to wrist (d)
2	16.2 – 18.2cm	14.7 – 16.2cm	10.9 – 11.7cm	24cm
4	18.2 – 20.2cm	16.2 – 17.7cm	11.7 – 12.5cm	27cm
6	20.2 – 22.2cm	17.7 – 19.2cm	12.5 – 13.3cm	30.5cm
8	22.2 – 24.2cm	19.2 – 20.7cm	13.3 – 14cm	34.0cm
10	24.2 – 26.2cm	20.7 – 22.2cm	14 – 14.7cm	38.0cm
12	26.2 – 28.2cm	22.2 – 23.6cm	14.7 – 15.4cm	42.0cm


COLOURS

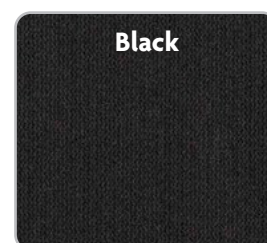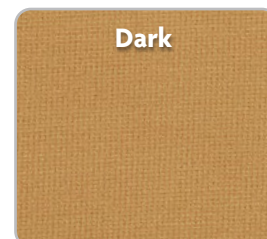


MEASUREMENTS

Measuring Tips:

- All circumferences and lengths to be measured to the nearest millimetre.
- Tape should be snug NOT tight.
- Arm should be placed at rest by side of body.

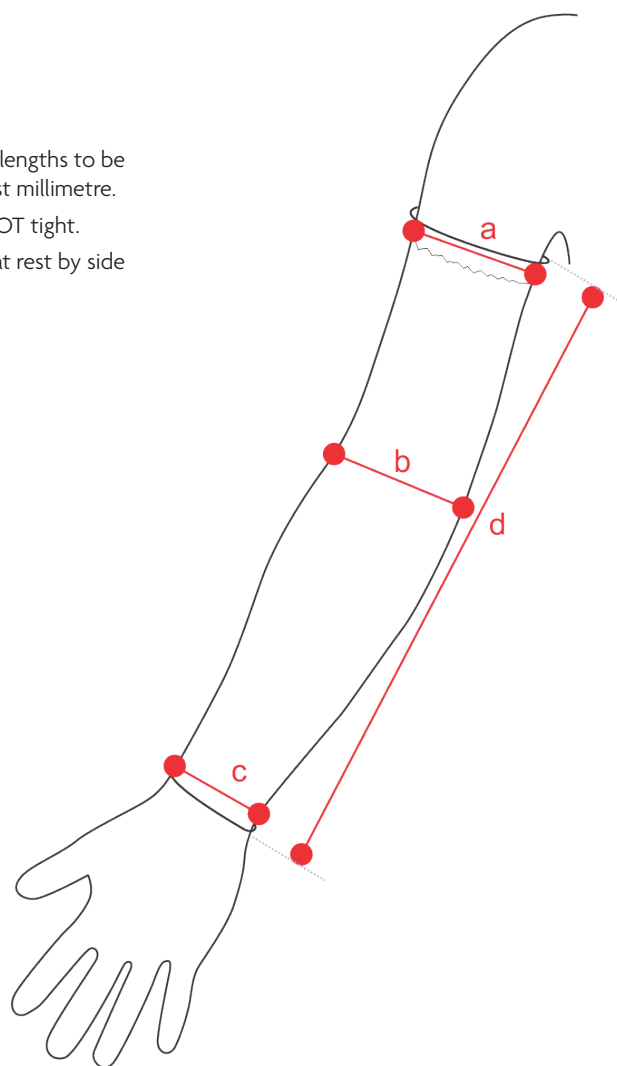

FABRIC OPTIONS

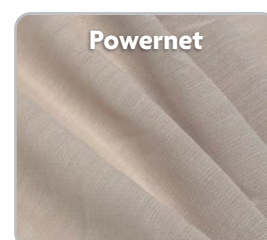

Please note that garments *cannot* be altered or repaired by Second Skin.



Second Skin Pty Ltd
40 O'Malley Street
Osborne Park WA Australia 6017
Email: perth@secondskin.com.au
Phone: +61 8 9201 9455



RF PAED-SC-UK 2026-05



Second Skin

CUSTOMISED COMPRESSION GARMENTS & SPLINTS

Fabric options

All Second Skin garments are made from a range of medical grade stretch fabrics. Several highly durable medical grade fabrics are incorporated in our garment designs, to address the individual needs of each client.

POWERNET FABRICS

Powernet is an open weave fabric which is prescribed for healed skin and for areas where additional support and compression is required. This fabric is particularly durable in work, school or outdoor activities. Powernet is often used in burns scar management and lymphoedema.

Our Powernet fabric now incorporates patented cooling properties, HeiQ Smart Temp, which provides outstanding thermal comfort to clients by keeping them cooler and dryer while wearing their garment.

Powernet fabric is available in light tan, dark tan and black.

SHIMMER FABRICS

Shimmer is a closed weave fabric which provides a softer interface with the skin. Shimmer fabric contours well and is comfortable and effective over joint surfaces and areas of skin fragility such as the elbow, ankle and knee. Shimmer fabric is used for early burns and flap management, external fixators, vascular conditions and certain types of lymphoedema.

Shimmer fabric is available in light tan, dark tan and black.

HYDROPHOBIC FABRICS

Hydrophobic fabric is used as both a lining and a base fabric in compression garments. This fabric wicks moisture away from the surface of the skin to evaporate, protecting the skin surface and ensuring optimal comfort.

Our hydrophobic fabric has an antimicrobial component embedded within the matrix of the fabric. This ensures our products are breathable and highly effective at removing bacteria from the fabric surface, to provide the wearer with all year round comfort. Due to this innovation, our garments are frequently prescribed to assist with the management of repeated skin breakdown and fragility.

Hydrophobic fabric is available in light tan, dark tan and black.

PERTH OFFICE

40 O'Malley Street
OSBORNE PARK WA 6017

P: +61 8 9201 9455

E: perth@secondskin.com.au

SYDNEY OFFICE

Suite 1, 8 Northcote Street
ST LEONARDS NSW 2065

P: +61 2 9386 0812

E: sydney@secondskin.com.au

BRISBANE OFFICE

Suite 1/41 Manilla Street
EAST BRISBANE QLD 4169

P: +61 7 3084 4319

E: brisbane@secondskin.com.au

MELBOURNE OFFICE

Unit 9, 202 Ferntree Gully Road
NOTTING HILL VIC 3168

P: +61 3 8637 0297

E: melbourne@secondskin.com.au